

 **Fasenra Pen**[™]
(benralizumab) Subcutaneous
Injection 30 mg

**USING YOUR
FASENRA PEN[™]**



(not actual size)

Before using your FASENRA PEN™, your healthcare professional will show you or your caregiver how to use it correctly. Read the **Instructions for Use** leaflet included with your medicine before you start using FASENRA PEN™ and each time you get a refill.

PREPARE




A. Gather supplies

Take the FASENRA PEN™ out of the refrigerator and remove it from the carton. Check the expiry date on your FASENRA PEN™. Do not use it if the expiry date has passed. Do not remove the cap until you are ready to inject. Set it on a flat surface and let it warm up at room temperature for **30 minutes**.

While you wait, you can gather the things you need for your injection, such as:




Alcohol wipes

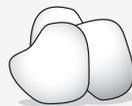

Cotton ball or gauze




Puncture-resistant sharps container

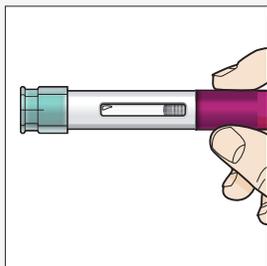

B. Check the medication

Check the liquid in your FASENRA PEN™ through the viewing window to make sure it's clear and colourless to yellow. It may contain small white particles. Don't use the FASENRA PEN™ if the liquid is cloudy, discoloured or contains large particles.


PREPARE

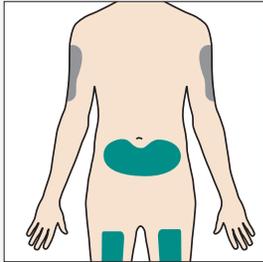

- Caregiver site only
- Self-injection or by caregiver

C. Choose the injection site


- You can choose one of three injection sites:
 - The front of your thigh (recommended)
 - Your lower abdomen
 - Your upper arm (only if a caregiver is injecting)
- Each time you inject, choose a different site that's at least 3 cm away from where you last injected. Do not inject within 5 cm of your belly button, where skin is tender, bruised, scaly or hard, into scars or damaged skin or through your clothing.
- Wash your hands with soap and water and clean the injection site using an alcohol wipe, letting the area air dry. Do not touch the cleaned area before injecting and do not fan or blow on the cleaned area.

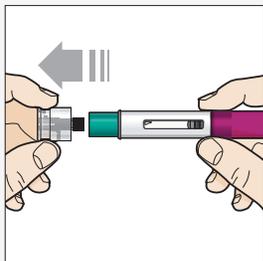


D. Pull off the cap

Carefully pull off the cap, and make sure you don't touch the needle or push on the green needle guard. Do not try to put the cap back on the FASENRA PEN™. You could cause the injection to happen too soon or damage the needle.


INJECT

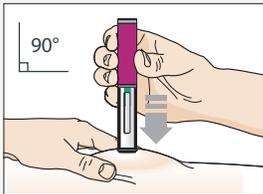

- You may gently pinch the skin at the injection site if you find that helpful.
- Line up the green needle guard at a 90-degree angle to the site you've chosen.

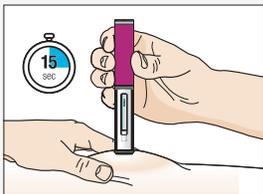

- Press down firmly. You will hear a **click**, which means that your injection has begun. The green plunger will start to fill the viewing window.
- Hold your FASENRA PEN™ in place for 15 seconds and do not change its position.
- You'll hear a **second click**, which means the injection is done. The green plunger will fill the viewing window.

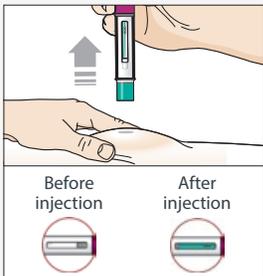

- Lift your FASENRA PEN™ straight up. The needle guard will slide down and lock into place.
- Check the viewing window to make sure you've injected all the liquid. If the green plunger doesn't fill the viewing window, you may not have received the full dose. If this happens, talk to your healthcare professional.
- If any bleeding occurs, use a cotton ball or gauze to gently hold pressure on the injection site until the bleeding stops. Do not rub the injection site.

Do not try to put the cap back on the FASENRA PEN™. Dispose of your used FASENRA PEN™ in a sharps container. Do not put a used pen or sharps container in your household waste or recycling.



Watch how to inject your FASENRA PEN™

FREQUENTLY ASKED QUESTIONS:

FASENRA PEN™ SELF-ADMINISTRATION

We have answers below to the most frequently asked questions about FASENRA PEN™. And remember, you can also talk to your doctor or pharmacist, or call 1-833-360-CONNECT (833-360-2666).

How to use FASENRA PEN™



How do I keep track of my injections?

Building a routine for taking FASENRA PEN™ can help keep track of your dosing schedule. The Connect360° Program can help you keep track of your injections by sending reminders by email or text.

You can also set up dose reminders and a [dose calendar](#) to help you stay on track with your treatment regimen.



Who do I contact if I need help with FASENRA PEN™ self-administration?

The Connect360° Program is here to answer any questions you have about your FASENRA PEN™. The Connect360° Program provides training for patients and caregivers injecting FASENRA PEN™. Speak to one of our dedicated nurses or schedule an appointment by calling 1-833-360-CONNECT (833-360-2666), Mon–Fri, 8 am–8 pm EST.

You can also watch the [Step-by-Step Video](#) on how to use your FASENRA PEN™.



Can someone else give me my injection?

Yes, a **trained** caregiver or healthcare professional can give your injection if you'd like. You and your doctor can make a decision about it together.

The Connect360° Program can provide training for patients and caregivers injecting FASENRA PEN™.



What if I miss a dose?

If you accidentally miss a dose of FASENRA PEN™, contact your healthcare professional and contact the Connect360° Program at 1-833-360-CONNECT (833-360-2666), Mon–Fri, 8 am–8 pm EST. They can advise you on when to take your next dose.

If you're worried about missing a dose, set up dose reminders and a [dose calendar](#) to help you stay on track with your treatment regimen.

FASENRA PEN™ packaging



What should I do if the expiration date has passed or if the medicine looks cloudy?

Don't use an expired FASENRA PEN™ or a pen where the liquid is cloudy, discoloured or contains large particles. Call the Connect360° Program at 1-833-360-CONNECT (833-360-2666), Mon–Fri, 8 am–8 pm EST, for further instructions.

Note that a small air bubble in the liquid is normal, and you can use your pen as usual.



What if my FASENRA PEN™ or carton looks damaged?

Do not use your FASENRA PEN™ if the security seal on the carton is broken or if the FASENRA PEN™ looks damaged. Call the Connect360° Program at 1-833-360-CONNECT (833-360-2666), Mon–Fri, 8 am–8 pm EST.

Storing FASENRA PEN™



How do I store my FASENRA PEN™?

Store your FASENRA PEN™ in a refrigerator at 2–8°C in the carton until you are ready to use it. Let it warm up at room temperature (20–25°C) for about 30 minutes before using it.



Can I store my FASENRA PEN™ at room temperature?

You should store your FASENRA PEN™ in the fridge, but it can be kept at room temperature up to 25°C for a maximum of 14 days.



What do I do if I left my FASENRA PEN™ out of the refrigerator for more than 14 days?

Discard your FASENRA PEN™ if it's been out of the fridge longer than 14 days and call the Connect360° Program at 1-833-360-CONNECT (833-360-2666), Mon–Fri, 8 am–8 EST, to get a replacement.



What if I need to travel with my FASENRA PEN™?

Always talk with your healthcare professional before travelling with prescription medication. You may need special documentation or storage containers. If you're leaving the country, it's also important to check the latest travel advice and advisories at travel.gc.ca.



Can I reuse my FASENRA PEN™ or share it with someone else?

No. Never reuse or share your FASENRA PEN™.

Disposing of FASENRA PEN™



How do I dispose of my used FASENRA PEN™?

Put your used FASENRA PEN™ in a puncture-resistant sharps container right away after use. Don't throw a used pen or sharps container in your household waste or recycling. Dispose of the sharps container as instructed by your healthcare professional or pharmacist.



Where can I get my injection supplies?

The Connect360° Program can provide you with alcohol swabs, as well as a puncture-resistant sharps container so you can properly dispose of your used FASENRA PEN™.

Where can I get a sharps container?

To get a new sharps container, call the Connect360° Program at 1-833-360-CONNECT (833-360-2666), Mon–Fri, 8 am–8 pm EST, or contact your pharmacy. Once your sharps container is full, take it to your local pharmacy.

Looking for additional resources to help you on
your treatment journey with FASENRA®?

Visit fasenra.ca

For complete information on FASENRA PEN™, refer to the **Patient Medication Information** and **Instructions for Use** leaflets included with your medicine. If you have any questions on how to use FASENRA PEN™, talk to your doctor or pharmacist or call 1-833-360-CONNECT (833-360-2666).

FASENRA® and the AstraZeneca logo are registered trademarks of, and FASENRA PEN™ is a trademark of, AstraZeneca AB, used under license by AstraZeneca Canada Inc.
© AstraZeneca Canada Inc. 2024
CA-7892E

